




Dheepika Lab

User guide

How to use Dheepika Lab, step by step, for patients and physiotherapists. Includes a map of the whole app, the care loop, a camera set-up guide and answers to common problems.

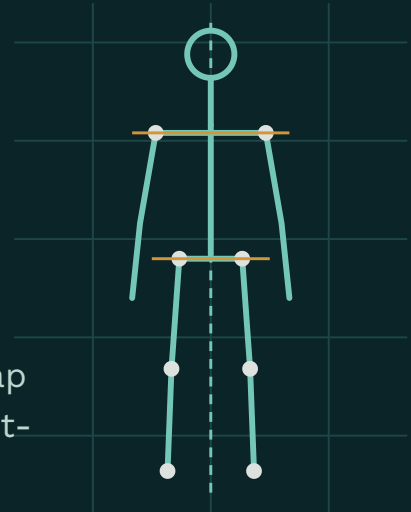

physiovision-ai-eta.vercel.app

For patients

Tell your physiotherapist where it hurts, do the movement tests with your phone camera, follow your exercise program and track your progress.

For physiotherapists

See every patient of the clinic, review camera estimates, sign the clinical impression and report, and build exercise programs.

Pilot software, not yet approved for patient care. Camera readings are estimates. They are not a diagnosis. Use test accounts only until clinical approval, patient consent, a data-processing agreement and real-device validation are in place. Screens in this guide show the built-in **simulated demo data**. Guide version 6 October 2026.

Inside this guide

- | | |
|--|---|
| <p>01 The app on one page
Map of every screen and how they connect</p> <hr/> <p>02 The care loop
How patient and physiotherapist work together</p> <hr/> <p>03 Getting started
Accounts, sign-in, demo mode</p> <hr/> <p>04 Patient: first-time setup
About you, consent, reason for visit</p> <hr/> <p>05 Patient: the Home screen
Today's session, check-in, navigation</p> <hr/> <p>06 Patient: doing an assessment
The six steps</p> <hr/> <p>07 Camera set-up guide
Where to put the phone and how to stand</p> <hr/> <p>08 Posture scan tools
Modes, views, self-correction, face cover</p> <hr/> <p>09 Patient: training, progress, profile
Exercise sessions and your data</p> <hr/> <p>10 Physiotherapist: the daily map
From new patient to reassessment</p> <hr/> | <p>11 Physiotherapist: overview and patients
Dashboard, directory, patient record</p> <hr/> <p>12 Reviewing an assessment
The seven workspace tabs</p> <hr/> <p>13 AI draft, "Why?" and the report
What is preliminary and what is signed</p> <hr/> <p>14 Programs, library and evidence
Building an exercise plan</p> <hr/> <p>15 Clinic: schedule and billing
Visits, attendance, courses, payments</p> <hr/> <p>16 Discharge summary
Closing an episode of care</p> <hr/> <p>17 Analytics and settings</p> <hr/> <p>18 How your data moves
Who sees what, privacy, deletion</p> <hr/> <p>19 Safety and limits
What the camera can and cannot do</p> <hr/> <p>20 Troubleshooting</p> <hr/> <p>21 Quick reference cards</p> <hr/> |
|--|---|

What Dheepika Lab is

It does


- Record where it hurts and the patient's own answers, exactly as given.
- Ask safety questions and stop routine tests if a warning sign is reported.
- Use the phone or laptop camera to estimate 33 body points, joint angles and posture lines. This runs **on the device**.
- Count exercise repetitions and show the target range.
- Give the physiotherapist a structured review, a report and a program builder.

It does not

- Diagnose any condition or name the injured tissue.
- See bones, joints or vertebrae inside the body. The camera sees only the body surface; there is no X-ray view.
- Measure centimetres from a single camera.
- Store video. Frames are processed in memory and discarded.
- Replace the physiotherapist. Every estimate waits for clinician review.


The app on one page

Everyone starts at the same website. The clinic decides who is a physiotherapist; everyone else gets the patient app. Both apps share records through the clinic server.

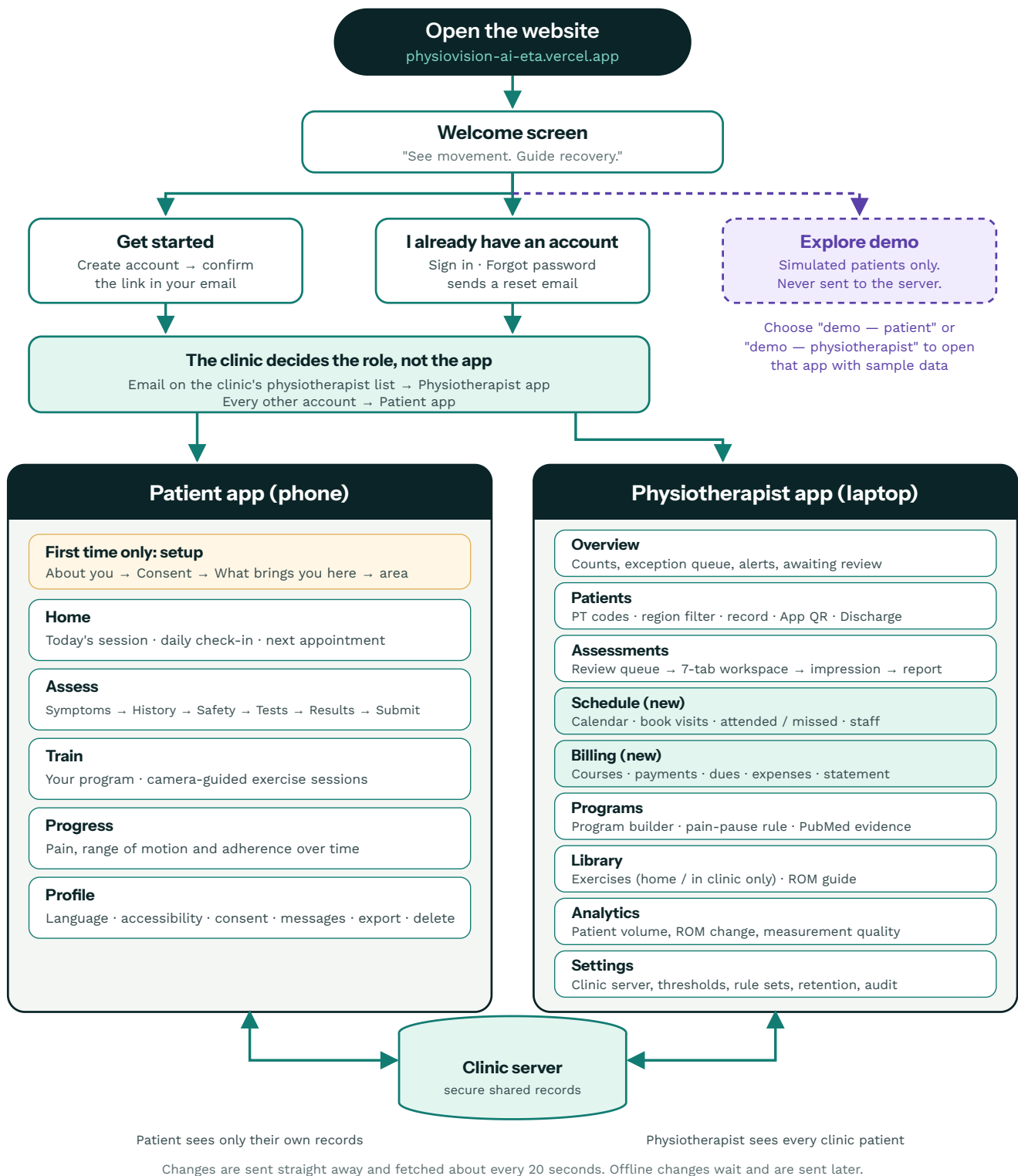

The care loop

Dheepika Lab follows one repeating loop. The patient provides answers and camera captures; the physiotherapist reviews them and makes every clinical decision.

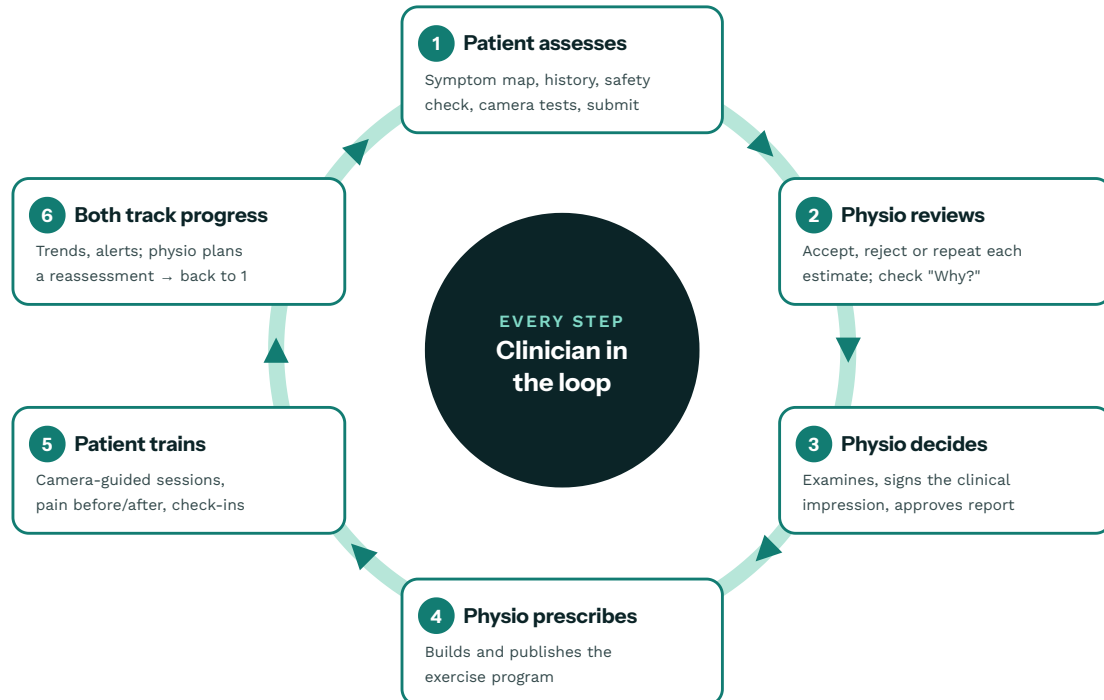

What the patient does

Answers questions honestly, does the camera tests in a safe space, trains as prescribed, and checks in. Stops if pain is sharp or worsening.

What the physiotherapist does

Reviews every estimate, examines the patient, writes and signs the impression, approves the report and sets the program.

What the app does

Checks capture quality, withholds numbers it cannot measure, flags threshold crossings for review, and keeps a full audit log.

Four kinds of information, always kept apart

Throughout the app, each piece of information carries a label that says where it came from. Look for these tags:

TAG	MEANING
Patient-reported	The patient's own answers: pain scores, symptom areas, history. Never changed by the app.
Camera-estimated	A number calculated from body points seen by the camera, only when the capture passed the quality checks.
Movement observation	An algorithmic observation : a camera estimate crossed a threshold set by the clinic. It is a flag for review, not a diagnosis.
Clinician-measured	Measured or decided by the physiotherapist, for example with a goniometer, or a signed impression.

Accounts and signing in

Create a patient account

- 1 Open **physiovision-ai-eta.vercel.app** in Chrome or Safari on your phone.
- 2 Tap **Get started**. Enter your name, email and a password, then tap **Create account**.
- 3 Open the email from the clinic and tap the confirmation link. It brings you back to the app, signed in.
- 4 Complete the first-time setup (next section).

Link expired or already used? The sign-in screen says so. Go back to **I already have an account** and sign in, or request a new link.

Sign in and passwords

- **I already have an account** → email and password → **Sign in**.
- **Forgot password?** sends a reset email. The link opens a screen to choose a new password.
- **Sign out** is at the bottom of the menu. It removes this account's records from the device, which keeps a shared phone private. Unsent changes wait for your next sign-in.

Physiotherapist accounts


You cannot choose "physiotherapist" in the app. The clinic owner adds the physiotherapist's email to the clinic's list on the server. That person then signs up or signs in with that email and gets the physiotherapist app. An account created earlier is upgraded automatically.

Demo mode

Explore demo — patient and **Explore demo — physiotherapist** open the app with **SIMULATED** patients. A purple banner shows on every screen. Demo data stays on the device and is never sent to the clinic server. Use demo mode to learn the app or to train staff.

Tip: add it to your home screen

In the browser menu, choose **Add to Home screen** (Android) or **Share → Add to Home Screen** (iPhone). The app then opens like any other app.

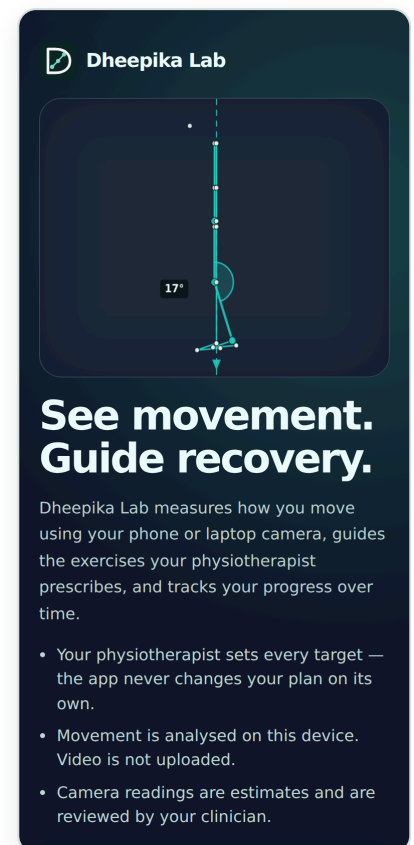


Welcome screen

Create account

First-time setup

A short setup runs once, right after your first sign-in. The progress bar shows "Step n of total".

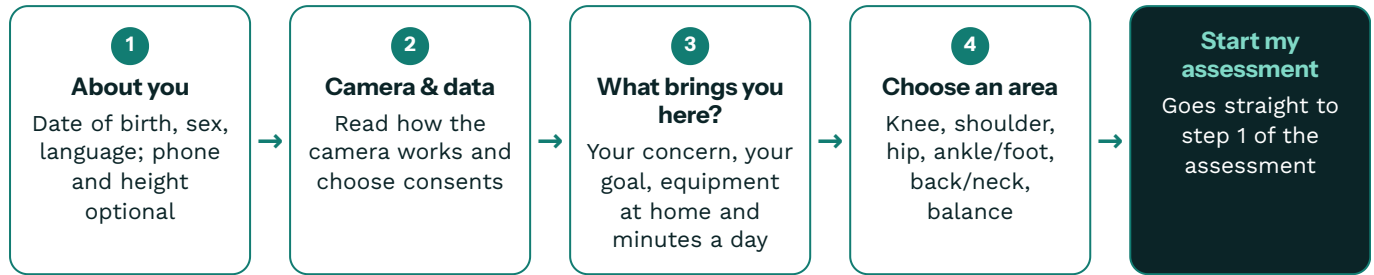

Step 3 also asks what equipment you have at home and how many minutes a day you can give; your physiotherapist sees both.

Your consent choices

The first two are needed for movement assessment. The last two are optional. You can change every choice later in **Profile → Privacy & consent**.

CHOICE	WHAT IT MEANS	NEEDED?
On-device camera processing	The camera finds 33 body points and calculates joint angles on your phone. Video frames are processed in memory and thrown away.	Required
Store and share results	Your measurements and questionnaire answers are saved and shared with your physiotherapist. Every access and change is recorded in an audit log.	Required
Store a still image with each posture scan	Keeps one photo per posture view so your physiotherapist can check it. You can ask for your face to be covered.	Optional
Anonymised research use	Anonymised measurements may be used to test and improve accuracy.	Optional


No video is stored. Only body-point coordinates and calculated measurements are saved, plus stills if you choose the optional image consent.

Which area first?

AREA	CAMERA TESTS AVAILABLE
Knee	Supported knee flexion (heel slide) · Five-times sit-to-stand · Double-leg squat (front view)
Shoulder	Active shoulder flexion (side view) · Active shoulder abduction (front view)
Hip	Active hip flexion, standing (side view) · Active hip abduction, standing (front view)
Ankle or foot	Knee-to-wall lunge (side view) · Double-leg heel raise (side view)
Back or neck	Standing forward bend (side view) · Standing side bend (front view) · Neck forward and backward movement (side view)
Balance and walking	Single-leg stance, eyes open (front view) · Marching on the spot (front view)

Every area also offers the optional standing posture scan. Your physiotherapist may change which tests you do. Later, use **Home → Assess another area**.

The Home screen

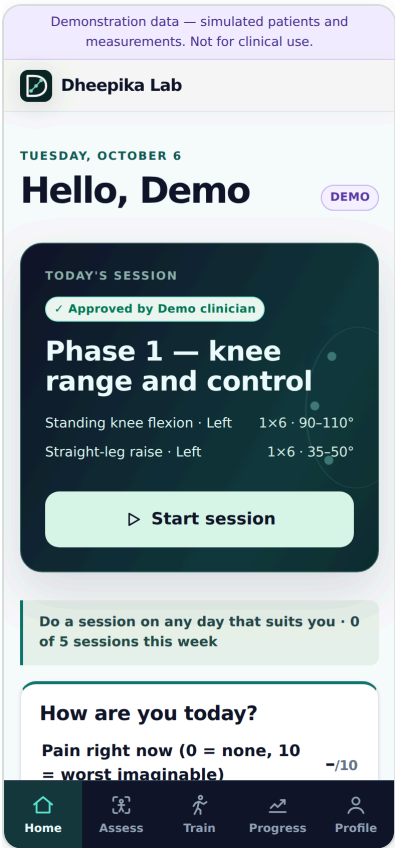

Home, with a program approved

What you see

- **Today's session:** your approved program, with each exercise, the side, sets × reps and the target range. Tap **Start session**.
- **This week:** sessions done out of the planned number. If you missed a day, just carry on with the next planned session; don't do extra to catch up.
- **How are you today?:** a daily check-in with your pain now (0–10) and an optional note. It goes to your physiotherapist. The app does not check it, so contact your clinic or urgent care if you feel unwell.
- **Next appointment, Assessment** status and **Recent sessions**.


No program yet? Home says so and asks you to complete an assessment, so your physiotherapist can build your plan.

The bottom menu

TAB	USE IT TO
Home	See today's plan, check in, and see appointments
Assess	Start, resume or view an assessment
Train	See your full program and session history
Progress	View pain, range-of-motion and adherence charts
Profile	Change language, accessibility, consent, messages, export or delete

Doing an assessment

Six steps. You can stop at any time and tap **Resume assessment** later. The bar at the top shows where you are (for example "2/6").

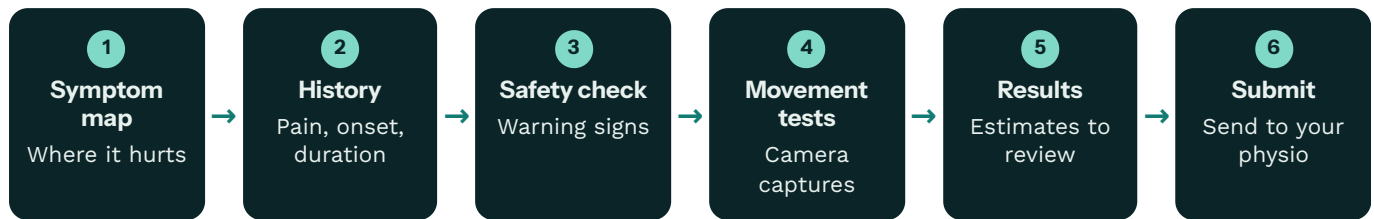


Step 1· Symptom map

- Turn the body with **Front / Right side / Back / Left side**. Labels show **your** left and right.
- Tap every area where you feel something, then choose what you feel there: pain, stiffness, weakness, numbness or tingling.
- For the knee, choose the exact spot: front / kneecap, inner side, outer side, back of knee or whole knee.
- If a symptom spreads, **draw its path** on the body.
- Prefer a list? Use **Choose areas from a list**.

The drawn body is an illustration to locate symptoms. It is not a camera image or a diagnosis.

Step 2· History

- Pain right now and worst in the last 24 hours (0–10).
- What it feels like, how long you have had it, how it started, what makes it worse, and whether it is constant.
- Questions adapt to your answers. They are saved **exactly as you give them**.

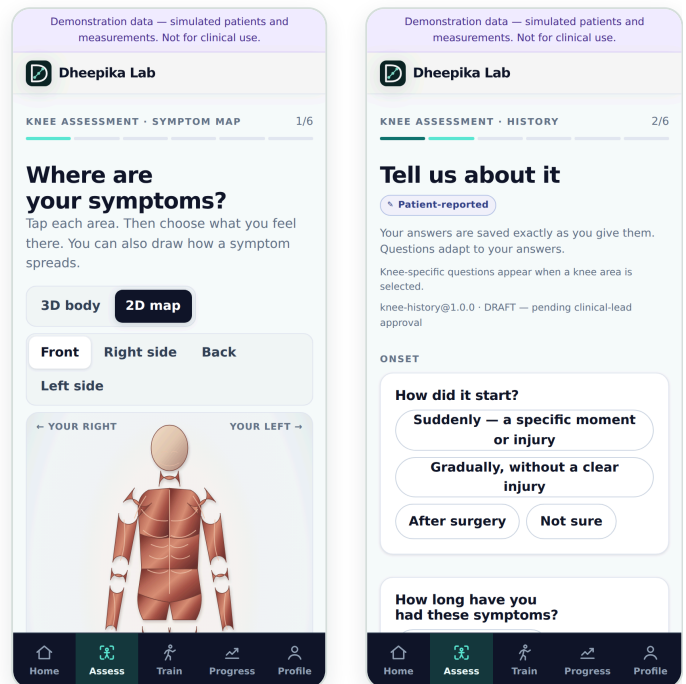

1 · Symptom map

2 · History

Step 3 · Safety check

Answer each question honestly. These questions help make sure exercise is appropriate right now. The answer decides what happens next:


RESULT	WHAT HAPPENS
Seek urgent care	"Please seek urgent medical care now." Tests stop. Your physiotherapist sees a critical alert.
Speak with your physio first	Tests are paused until your physiotherapist reviews your safety answers.
Clear	Continue to the movement tests.

In an emergency, do not wait for the app. Call your local emergency number. Profile shows the clinic's emergency number.


Step 4 · Movement tests

Your test plan is listed with a **Start** button for each test. Each test explains its own camera set-up and counts down before recording. A test shows ✓ **done** when the capture was valid, or ✗ **redo** when it was not. You can **Recapture** any test. The optional **static standing posture** scan (front, side and back) is at the bottom of the list.


Before your first camera test, read the **Camera set-up guide** on the next page.

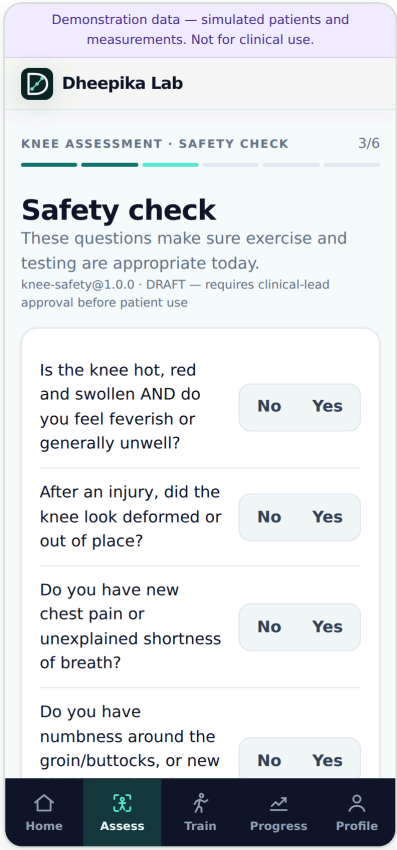


3 · Safety check

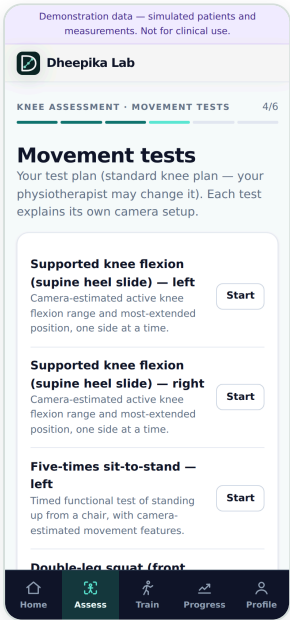

4 · Movement tests

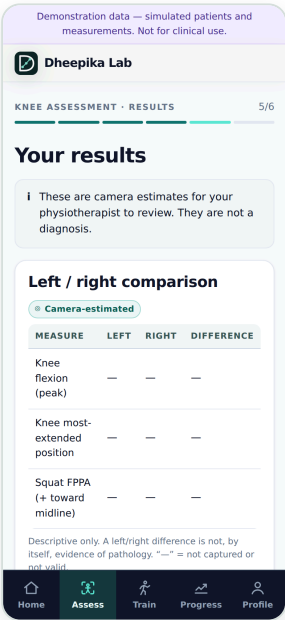

5 · Results

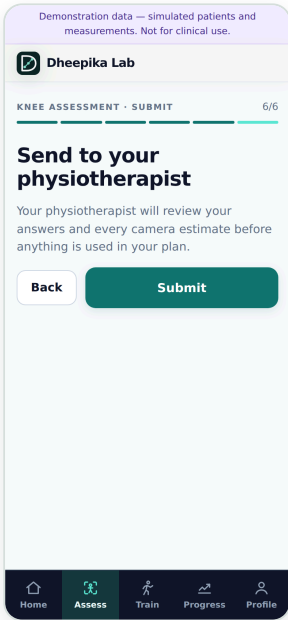

6 · Submit

Step 5 · Results


Your camera estimates, such as a left/right comparison. A dash "—" means "not captured or not valid". These are estimates for your physiotherapist to review. They are not a diagnosis.

Step 6 · Submit

Tap **Submit to my physiotherapist**. Home then shows "Submitted — awaiting physiotherapist review", then "Reviewed by your physiotherapist" when they have checked it.

Camera set-up guide

Good set-up is the difference between a valid measurement and a "—". Each test tells you where to put the phone: usually 2–3 m away (1.5–2 m for ankle tests), at knee, hip or chest height, facing you or to your side.

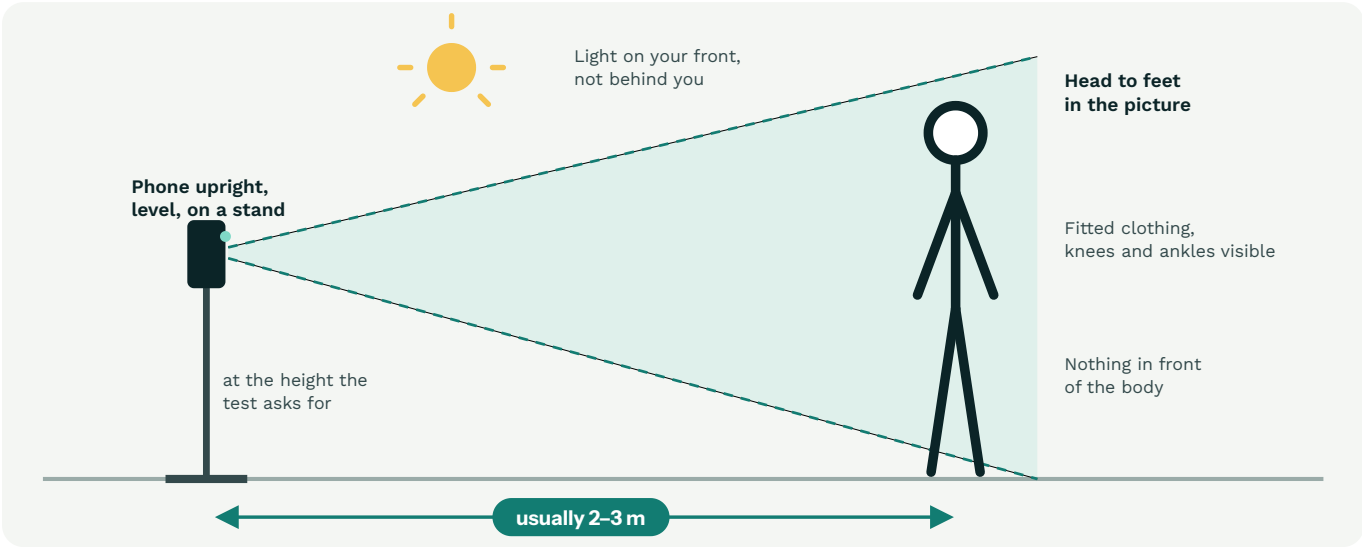


The set-up checklist

Before recording, the app checks the set-up. Each line turns green when it passes. The capture waits until all of them pass.

CHECK	IF IT DOES NOT PASS
Body detected	Step into view. Only one person should be in the picture.
Required body parts in frame	Make sure the joints this test needs are inside the picture and not covered.
Distance	The line shows how much of the picture height your body fills and the range to aim for, for example <i>"body fills 78% (aim for 30–97%)"</i> . Move closer or further back to get into the range.
Body orientation	Face the direction the test asks for: front, left side or right side.
Camera level	Stand the phone upright and straight, not tilted.
Lighting	Turn on a light in front of you. Avoid a bright window behind you.

The checks are smoothed, so a small sway does not flip "ready" on and off. The distance range adapts to your camera's resolution. On a sharp HD phone camera you can stand further back.

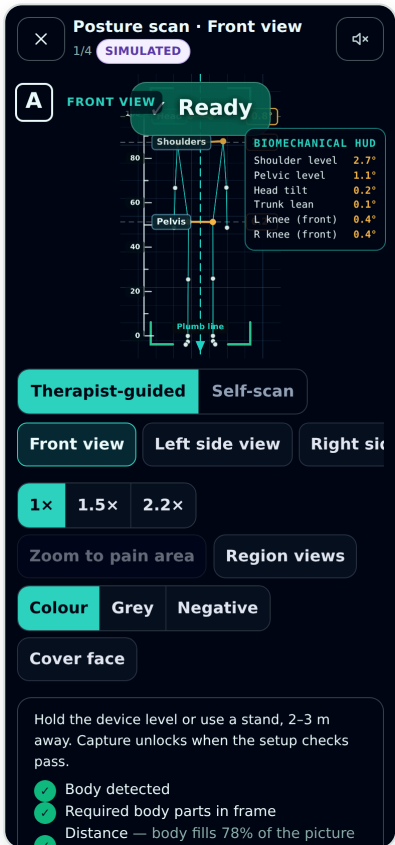

Posture scan tools

The standing posture scan takes front, side and back views, with a spoken instruction for each. Each view is a 3-second capture; the app uses the middle (median) value.

Two scan modes

Therapist-guided: the physiotherapist holds the device or uses a stand, 2–3 m away, and uses the rear camera. **Capture this view** unlocks only when the set-up checks pass.

Self-scan: the patient places the phone on a stand. The capture starts automatically when the checks pass and you hold still.

What the lines mean

- **Ruler:** percentage of ankle-to-nose height (0 = ankles, 100 = nose). It is a relative scale, not centimetres.
- **Plumb line:** dashed vertical line from the ankles.
- **Level lines:** head, shoulders and pelvis, with the measured angle.
- **A / L / R / P badge:** which view is being captured.

Display tools

TOOL	WHAT IT DOES
1x · 1.5x · 2.2x	Zooms the display only. The measurement is the same at every zoom.
Zoom to pain area	Zooms to the area the patient marked on the symptom map and outlines it as "Patient-reported". The button is off if no area was marked, or if that part of the body is not in view.
Region views	Shows four zoomed views at once (full body, head & neck, trunk, lower limb). They are crops of the same single camera, not extra cameras.
Colour · Grey · Negative	Changes the image tone for easier viewing. It shows the body surface only; it is not an X-ray .
Cover face	Hides the face on screen, and on saved stills.

Why a value shows "—"

A number appears only when every body point it needs is clearly visible. Otherwise the app withholds it and says why ("Waiting" means the set-up checks are not finished yet):

MESSAGE	WHAT TO DO
Hands or an object in front of the body	Arms relaxed by your sides; put the phone or any object away.
... not in view	Step back or re-aim the camera so the named joints are inside the picture.
... not clearly visible / covered	Loose clothing or shadows hide the joint. Wear fitted clothing and add light.
Not steady or clear enough during the capture	Stand still for the full 3 seconds and try again.
Could not be calculated in this position	Check you face the right way for this view.

A partly hidden body never gives a full-body result. If the phone or an arm covers part of the body, only the measurements that can be seen are shown.

Self-correction: before and after


After a view is captured, you can tap **Try a 15-second self-correction**. A countdown starts and the app speaks the instruction: stand as tall and balanced as feels comfortable, let your shoulders settle, then hold still.

The view is then captured again and a table shows **Before, After** and **Change** for each measure.

- Only measures taken in **both** captures get a number. Anything else shows "—".
- This is for coaching. Your usual posture (the first capture) is the one saved for your physiotherapist.
- Tap **Self-correct again** to repeat, or **Retake** to start the view over.


Spoken prompts

With the voice coach on, the app tells you how to stand for each view: face the camera; left side to the camera; right side to the camera; back to the camera.

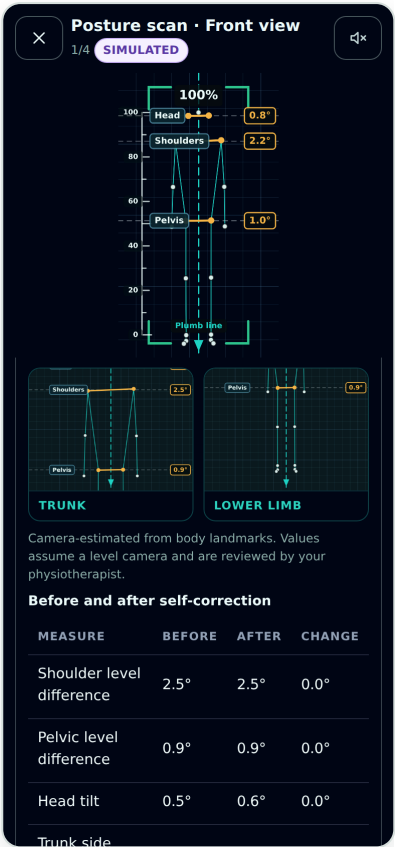


Before / after table

Training, progress and profile

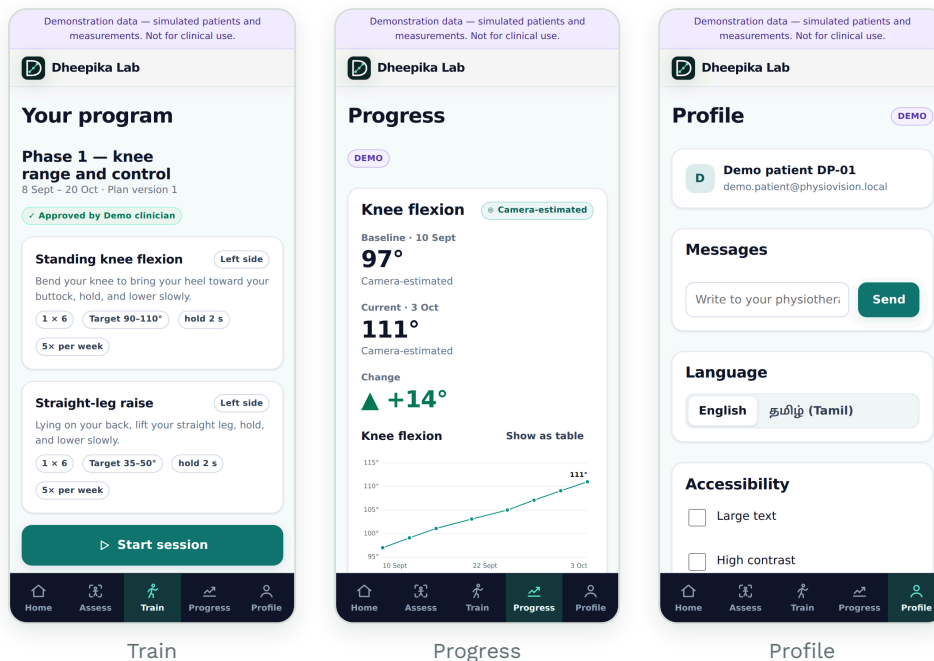

An exercise session

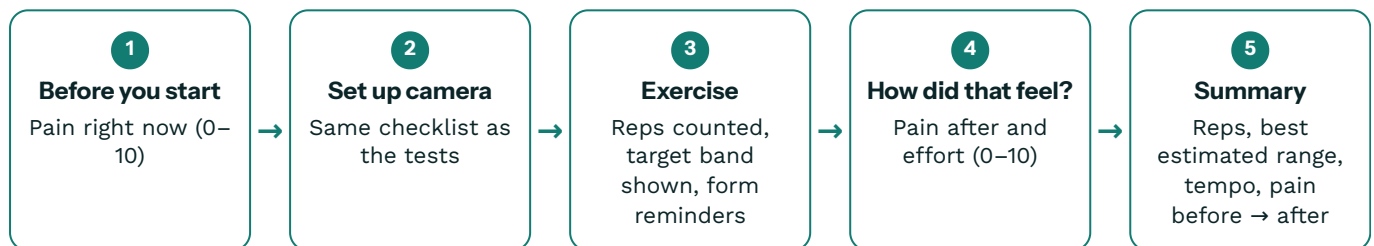

- The **Train** tab shows the program, who approved it, sets × reps, target range, hold time, how often per week, and your session history. Exercises that do not need the camera are listed under **Also in your plan**.
- If your pain rises past the limit your physiotherapist set, the session pauses and your physiotherapist is told. **Stop if pain is sharp or getting worse.**
- If a session is interrupted, the progress so far is saved.

Progress


Charts of pain (0–10), range of motion and adherence, from baseline to now. **Camera-estimated** values (session peaks) and **clinician-measured** values (goniometer) are drawn differently, so you can tell them apart. Tap **Show as table** for exact numbers.

Profile

- Language (English; Tamil is a draft awaiting review).
- Large text, high contrast, reduced motion, voice coach, captions, vibration.
- Privacy & consent; the clinic's emergency number; messages to your physiotherapist.
- **Export my data**, or **Delete my account and data** (type DELETE to confirm). Deleting removes your records from the clinic server and the device.

The physiotherapist's map

A typical path from a new patient to a reassessment. Use a laptop or tablet for review work. On a phone the menu shows Overview, Patients, Assessments, Schedule and **More** for everything else.

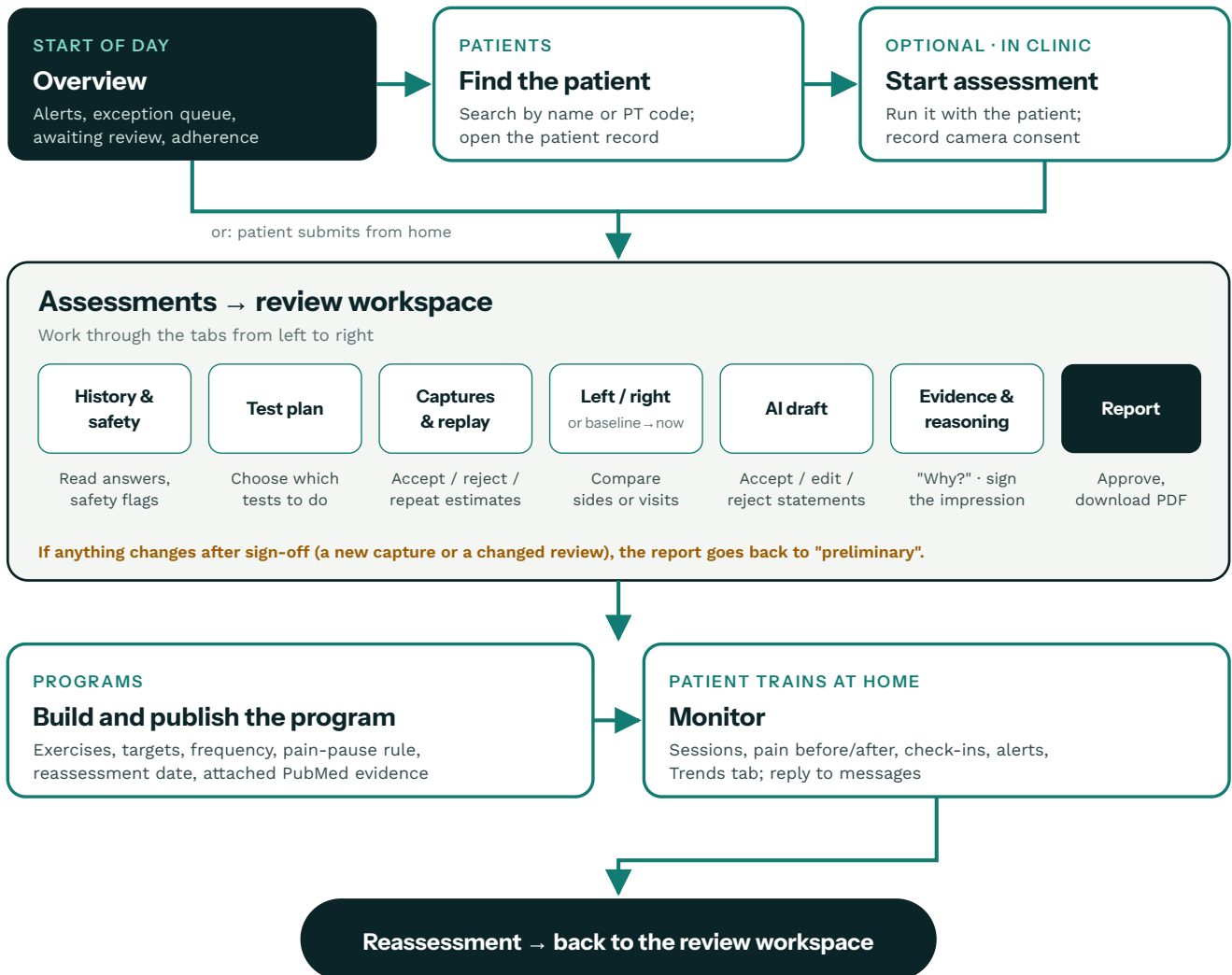

Overview and patients

Dheepika Lab

Overview

Patients

Assessments

Schedule

Billing

Programs

Library

Analytics

Settings

Sign out

Demonstration data — simulated patients and measurements. Not for clinical use.

DHEEPIKA LAB

Overview

Demo clinician · Physiotherapist (simulation)

+ Programs

Active patients

1

5 total

Assessments awaiting review

5

Requires clinician review >

Sessions (7 days)

2

Open alerts

5

1 critical

Exception queue

Rules (6 unreviewed)

Demo patient DP-01

4 of 10 planned sessions in the last 14 days

plan Phase 1 — knee range and control v1, 5/week

rule threshold unreviewed

Reported symptom changes & alerts

Demo patient DP-03

Red flag — urgent · urgent: 6 Oct

calf (knee-safety@1.0.0)

Awaiting review [View all](#)

Demo patient DP-01

4 Oct, 02:28 pm · 11

pending estimates

Program adherence (28 days)

Demo patient DP-01

40%

low

Overview (demo data)

Overview

- **Counts:** active patients, assessments awaiting review, sessions in the last 7 days, and open alerts (critical ones counted separately).
- **Exception queue:** patients who need attention, such as low adherence against the plan.
- **Reported symptom changes & alerts:** safety flags, such as a red flag that needs urgent care.
- **Awaiting review** and **Program adherence (28 days)**.

Patients directory

Every patient of the clinic, from every device, with a **PT code** (for example PT-3F9A21), age and sex, main complaint and registration date. Search by name or PT code, and filter by body region (knee, shoulder, hip, ankle, back and neck, balance). Use **Start assessment** to run an assessment with the patient in the clinic. You pick the area, and if the patient has not yet given camera consent you can record it there.

Overview

Patients

Assessments

Schedule

Billing

Programs

Library

Analytics

Settings

Sign out

Demonstration data — simulated patients and measurements. Not for clinical use.

Registered patient directory

5 patients

Search by name, patient ID or coi

All body regions

AllNeeds reviewAlertsLow adherence

PATIENT ID	PATIENT · AGE / SEX	PRESENTING COMPLAINT (PATIENT'S WORDS)	REGISTERED	PROGRAM	ADHERENCE	ALERT
PT-967330	DP Demo patient DP-01 34 y · male	-	27 Aug 2026	Phase 1 — knee range and control last session 3 Oct	40%	! 1
PT-DD62BD	DP Demo patient DP-02 58 y · female	-	27 Aug 2026	None	-	! 1
PT-FA4AEC	DP Demo patient DP-03 45 y · female	-	27 Aug 2026	None	-	! 1
PT-891E9E	DP Demo patient DP-04 51 y · female	-	27 Aug 2026	None	-	! 1
PT-E9785B	DP Demo patient DP-05	-	27 Aug 2026	None	-	! 1

Patients directory

The patient record

Header buttons: **Start assessment**, **Build program**, **Discharge** and **App QR**. The QR code opens the app's sign-in page on the patient's own phone; it contains only the app address. The header also shows the patient's equipment at home and daily time, as they reported them.

TAB	CONTENTS
Overview	Summary, current program, alerts
Assessment	Assessments for this patient and their status
Movement measurements	Camera estimates and clinician measurements, each with its source
Patient-reported outcomes	Questionnaire answers and pain scores
Programs	Current and past programs, with versions
Trends	Changes over time
Devices & export	Shared step and walking data, and data export
Session history · Progress	Every exercise session; progress charts
Clinical notes	Your notes on this patient

Reviewing an assessment

Open **Assessments**, pick a submitted assessment, and work through the seven tabs.

Overview

Patients

Assessments

Schedule

Billing

Programs

Library

Analytics

Settings

Sign out

Demonstration data — simulated patients and measurements. Not for clinical use.

KNEE ASSESSMENT · SUBMITTED

Demo patient DP-02

DEMO

58 y · started 5 Oct, 01:29 pm · submitted 5 Oct, 02:29 pm

History & safety

Test plan

Captures & replay

Left / right

AI draft

Evidence & reasoning

Report

Showing the latest capture per test (4 captures in total, earlier attempts retained). Accepting a camera estimate records that you reviewed it; it remains labelled camera-estimated.

Supported knee flexion (supine heel slide) — left

knee_supported_flexion@1.1.0 · pv-knee-1.1.0 · synthetic-skeleton 1.0.0 · view lateral_left · 5 Oct, 01:29 pm

Camera-estimated

SIMULATED

valid capture

Knee flexion (peak)

124.1°

per rep: 124.1, 123.9, 123.8

Knee extension (most-extended position)

2.9°

0° = straight. Higher = knee did not fully straighten. Hyperextension cannot be measured by this 2D method.

Quality: coverage 100% · mean landmark confidence 0.97 · 3 valid of 3 attempts · 30 fps

Show synchronized replay

Reference measurement for validation (0 recorded)

Clinician-measured

knee flexion peak

pending

Accept

Reject

Repeat

knee extension position

pending

Accept

Reject

Repeat

Supported knee flexion (supine heel slide) — right


Captures & replay: estimate, version, view, quality line and review buttons

TAB	WHAT YOU DO THERE
History & safety	Read the symptom map, radiation path and the original answers. Safety flags are shown at the top. A safety hold withholds routine automated conclusions.
Test plan	Choose which tests the patient should do. The patient's app follows your plan.
Captures & replay	For each test, see the latest valid capture with its protocol version, camera view, quality (coverage, confidence, valid attempts, frame rate) and a synchronised replay. Accept , Reject or Repeat each estimate. Enter a reference measurement (for example a goniometer reading) for validation. The posture board shows each view with region zooms, and a grey or negative tone.
Left / right	Compare sides. For a reassessment, this tab shows baseline → current. A difference is descriptive only, not evidence of pathology by itself.
AI draft	A structured draft (see next page) to accept, edit, reject or defer, statement by statement.
Evidence & reasoning	The evidence map and "Why?" view. Record your examination findings and the clinical impression .
Report	Preview, approve and download the report.

Dheepika Lab · User guide · Pilot, not for patient care

17 / 27

AI draft, "Why?" and the report

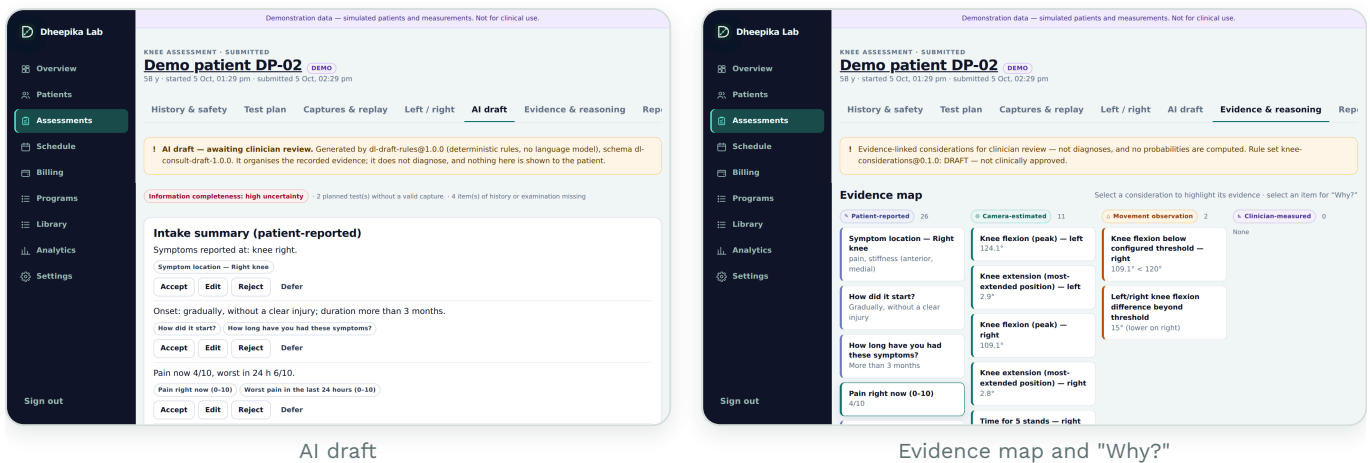


AI draft

The draft is built from fixed, versioned rules, not a language model. It organises the recorded evidence into five sections: intake summary (patient-reported), observations (camera-estimated / algorithmic), safety flags, contradictions to clarify, and considerations. Every statement shows the evidence it cites. For each one, choose **Accept**, **Edit**, **Reject** or **Defer**. The original proposal is kept next to your decision. An **information completeness** banner lists missing tests or history. Nothing in the draft is shown to the patient.

Evidence & reasoning: the "Why?" view

The evidence map has four columns: **Patient-reported**, **Camera-estimated**, **Movement observation**, and **Clinician-measured**. Select an item to see **why** it is there: the measurement, camera view, confidence, calculation, rule version and limitations. Considerations are evidence-linked prompts for your review. They are not diagnoses, and no probabilities are computed.

The report

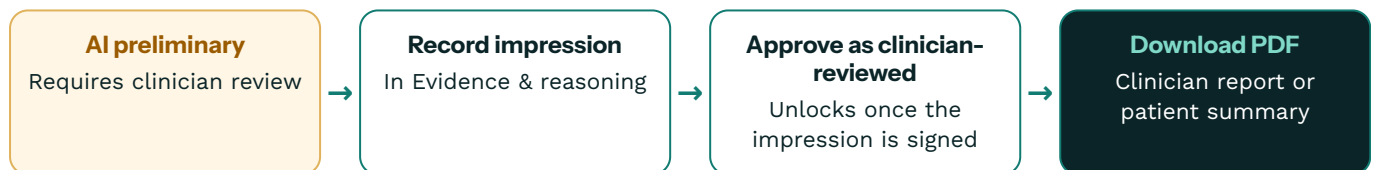

The report uses the latest valid capture for each test and view. If a capture or review changes after approval, the report goes back to **preliminary** until you approve it again. Use **Print view** to print.

Programs, library and evidence

Demonstration data — simulated patients and measurements. Not for clinical use.

Dheepika Lab

- Overview
- Patients
- Assessments
- Schedule
- Billing
- Programs**
- Library
- Analytics
- Settings
- Sign out

Program builder

Clinician-approved

Patient
Select...

Program title
Phase 1 — range of motion

Start date
10/06/2026

End date
11/17/2026

#1
Standing knee flexion (v1.0.0)
Left
Right
Remove

Bend your knee to bring your heel toward your buttock, hold, and lower slowly. Measures knee flexion (camera-estimated, lateral view). Supported target range 10–150°.

Sets	Repetitions	ROM target min (°)	ROM target max (°)
2	10	80	100

Hold (s)	Rest between sets (s)	Min. rep duration (s) — tempo	Frequency (per week)
3	45	2.5	5

Special instructions (shown to patient)

Program builder

Building a program

- 1 Open **Programs** → **new**. Choose the patient, program title, start and end dates.
- 2 Add camera-tracked exercises. For each one, set the side, sets, repetitions, range-of-motion target (min–max °), hold, rest, minimum rep duration (tempo), frequency per week and special instructions.
- 3 Set the **pain-pause rule**: pause the session if patient-reported pain reaches a level, or rises by a set number of points above the pre-session score. Leave it blank to switch it off.
- 4 Add **library exercises** that do not need the camera (sets, reps, hold or duration, frequency).
- 5 Set the **reassessment** date and pauses.
- 6 Optional: **Search PubMed** and **Attach to plan** supporting articles. Records are shown exactly as PubMed returns them; the app never writes citations itself.
- 7 Publish. The patient sees it on Home as "Approved by ...". Editing creates a new version with a list of changes from the previous one.

Library

The **Exercise library** holds exercises that do not need the camera. Each one is marked **Home after prescription** or **In clinic only**; in-clinic-only items cannot be added to a home plan, and changing the level creates a new version for review. The program builder flags items that need equipment the patient did not list. The **ROM measurement guide** (button on the Library page) shows goniometer placement for each joint the app measures, without normal ranges.

Schedule and billing

Run the clinic day in the app: book visits, mark who attended, and keep courses, payments and expenses in one place. Every figure is a total of what staff record; nothing is estimated.

Demonstration data — simulated patients and measurements. Not for clinical use.

Dheepika Lab

- Overview
- Patients
- Assessments
- Schedule**
- Billing
- Programs
- Library
- Analytics
- Settings

Sign out

CLINIC

Schedule

Book visits, record attendance and see each day's agenda. Session counts on treatment courses come from attendance recorded here.

Billing

October 2026 < Today >

MON	TUE	WED	THU	FRI	SAT	SUN
			1	2	3 1	4
5	6	7 1	8	9	10	11 2
12	13	14	15	16	17	18
19	20	21	22	23	24	25
26	27	28	29	30	31	

Saturday 3 October

Demo patient DP-01 PT-

10:00 BE24A4 **DEMO** **Attended** **Undo**

Treatment session · Demo clinician ·
Knee rehabilitation (demo) (3/8)

The WhatsApp reminder contains only the date, time, clinic and staff name. It opens WhatsApp on this device; nothing is sent until you press send there.

Book a visit

Patient	Date	Time	Staff	Type
Select... ▾	10/03/2026 📅	09:00 ▾	Demo clinician ▾	Treatment ses... ▾

Treatment course Note (optional)

Schedule

- 1 Pick a day on the calendar. The number under a date is how many visits are booked.
- 2 Under **Book a visit** choose the patient, date, time, staff member, visit type and, for treatment sessions, the patient's course. A second booking for the same staff member at the same time is blocked.
- 3 After the visit, mark it **Attended** or **Missed** (or **Cancel**). **Undo** reverses a mistake.
- 4 If the patient has a phone number, **WhatsApp reminder** opens WhatsApp with a message carrying only the date, time, clinic and staff name. Nothing is sent until you press send in WhatsApp.

Add or pause staff under **Staff**. Adding a name there does not create a login; physiotherapist logins are added by the clinic owner.

Overview

Patients

Assessments

Schedule

Billing

Programs

Library

Analytics

Settings

Sign out

CLINIC

Billing

Treatment courses, payments, expenses and monthly statements. Figures are totals of what staff record here.

Schedule

Received in Oct 26
₹0.00

Expenses in Oct 26
₹15,000.00

Net in Oct 26
-₹15,000.00

Outstanding dues
₹3,200.00
1 patient(s)

Treatment courses

Export CSV

Attended and missed are counted from visits marked on the schedule. Courses do not show recovery; that is judged by the physiotherapist from assessments.

PATIENT	COURSE	SESSIONS	AGREED FEE	STATUS
Demo patient DP-01 DEMO	Knee rehabilitation (demo) from 2026-09-22	3/8 attended · 1 missed · 2 booked	₹6,400.00	active

New course

Patient

Course name

Planned sessions

Agreed fee (₹)

Start date

Select...

e.g. Knee rehabilita

10/06/2026

Add course

Billing: this month, courses, dues, payments, expenses and statement

Billing

PART	WHAT IT DOES
Treatment courses	The sessions and fee agreed with the patient. Attended and missed counts come from the schedule. Mark a course Complete or Stop .
Outstanding dues	For each patient: agreed fees minus payments received.
Payments received	Amount in rupees, method (cash, UPI, card, bank transfer), date and an optional receipt or UPI reference.
Expenses	Rent, salaries, utilities, equipment, supplies, software and other costs.
Monthly statement	Money received, expenses and net for each of the last 12 months, as a chart and a table.
Export CSV	Courses, payments, expenses or the statement as a spreadsheet file. Each export is recorded in the audit log.

Courses do not show "recovery". Progress is judged by the physiotherapist from assessments. Only physiotherapist accounts can record courses, payments and expenses.

Discharge summary

Dheepika Lab

Overview

Patients

Assessments

Schedule

Billing

Programs

Library

Analytics

Settings

Sign out

Demonstration data — simulated patients and measurements. Not for clinical use.

DISCHARGE

Demo patient DP-01 DEMO

PT-BE24A4

Back to patient

Discharge summary Preliminary — not signed

Attendance, pain scores and measurements are filled in from the record. Write the summary in your own words; nothing here is generated.

Summary of treatment and outcome

Demo text: course completed; stairs reported easier by the patient.

Home advice

Follow-up

Sign discharge summary


Download PDF

Discharge summary: physiotherapist's text above, the assembled summary below

- 1 Open the patient and press **Discharge**.
- 2 The summary fills itself from the record: visits attended, the treatment course, the patient's first and latest pain scores, the latest assessment's measurements, and the signed clinical impression.
- 3 Write the **summary of treatment and outcome**, **home advice** and **follow-up** in your own words. Nothing in these boxes is generated.
- 4 Press **Sign discharge summary**, then **Download PDF**. Signed versions are kept and cannot be edited; signing again adds a new version.

It goes back to "preliminary" when data changes. If a visit is marked attended, a capture is added or a course changes after signing, the summary shows "Signed version out of date" until you sign again.

Analytics and settings

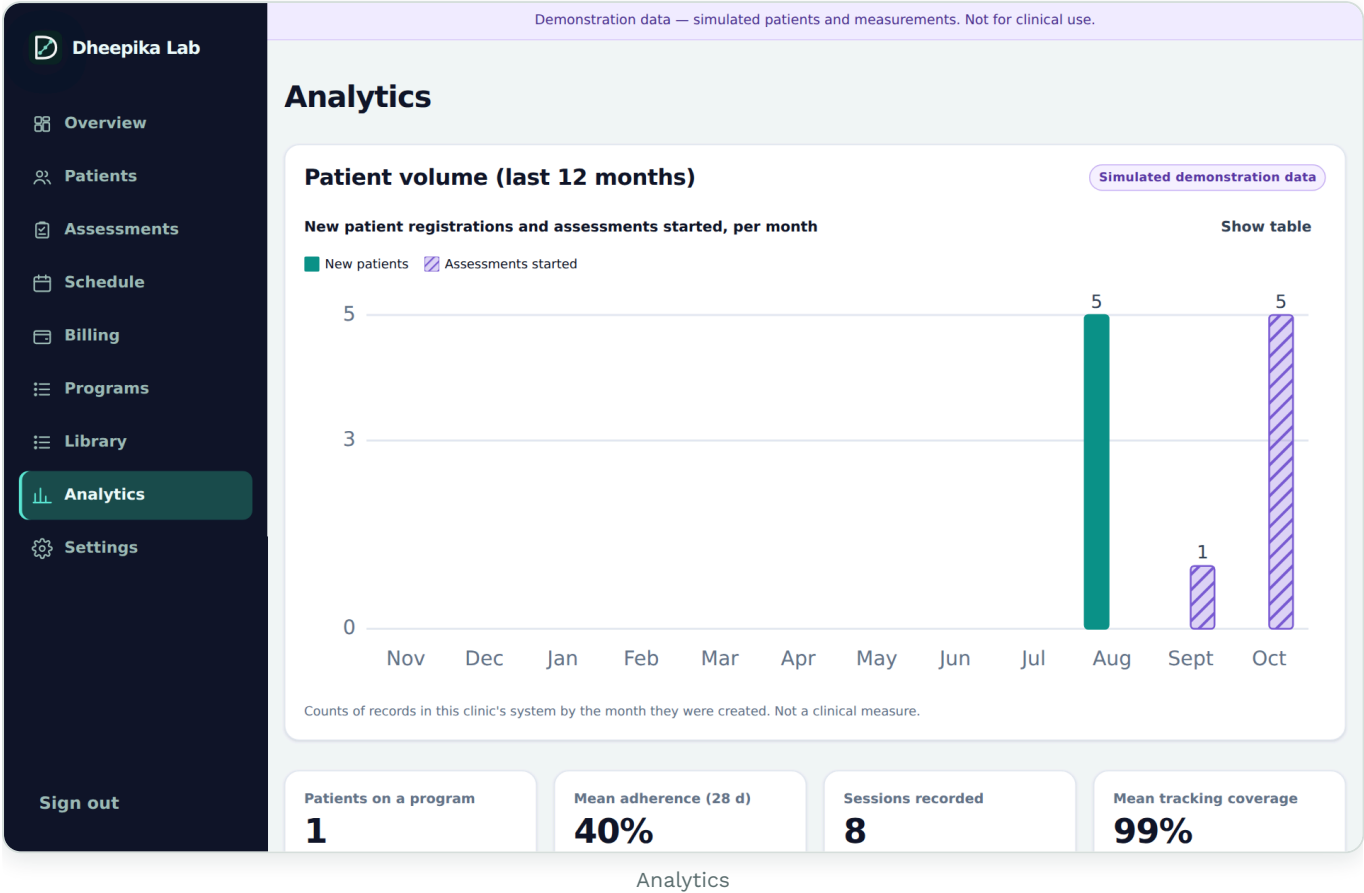

Analytics

- **Patient volume (last 12 months)**: a bar chart, with a table view.
- **Mean change in camera-estimated ROM** (first → latest session).
- **Measurement quality** of the captures.

Settings

SECTION	USE
Clinic server	Sync status ("Up to date", changes waiting, last synced) and Sync now .
Clinic	Clinic name, the emergency number shown to patients, and whether payment records are kept when a patient deletes their account (off by default). Release readiness and incidents on this device are listed below it.
Observation thresholds	The thresholds that turn a camera estimate into a movement observation. Unreviewed rules are counted on the Overview.
Clinical rule sets & protocols	Versions and approval status. Draft rule sets are labelled "not clinically approved".
Motion engine (this device)	Body-tracking model: Lite (fastest), Full , or Heavy (most accurate; a ~31 MB download, best on a laptop).
Demonstration data	How demo patients are flagged. Sign in with a real account to remove demo data.
Data boundary · Audit log	What stays on the device and what is shared; the latest 40 audit entries.

How your data moves

ACCOUNT	CAN SEE	CAN CHANGE
Physiotherapist	Every patient of the clinic, on every device	Everything in the clinic, except history, which can never be edited
Patient	Only their own records	Only their own records. Never plans, impressions, reviewed reports or clinic settings
Demo	Simulated data only	On the device only. Never sent to the server
Schedule and billing	Physiotherapists. Clinic expenses are never visible to patients	Physiotherapists only

What is saved

- Answers, pain scores and check-ins.
- Body-point coordinates and calculated measurements.
- Posture stills, only with the optional image consent.
- Reviews, impressions, programs and an audit log of every access and change.

What is never saved

- Video. Frames are processed in memory and discarded.
- Demo data on the server.
- Your password in the app. Sign-in is handled by the clinic server.

Syncing

- Every change is queued and sent to the clinic server. Changes from other devices are fetched about every 20 seconds, when the connection returns, and when you reopen the app.
- **Offline?** A banner says how many changes are waiting. They are kept on the device and sent automatically later.
- If the same record is changed in two places at once, the later change is kept.

Leaving

- **Sign out** removes this account's records from the device.
- **Profile → Delete my account and data** permanently removes the patient's records and sign-in from the server and the device. If the clinic has switched on **Keep payment records** (Settings → Clinic, off by default), payment amounts, dates, method and receipt reference are kept for accounting with the patient link, course name and notes removed.
- **Profile → Export my data** downloads a copy.

Safety and limits

Not for patient care yet. Dheepika Lab is a pilot. Clinical approval is on hold. Use test accounts and consenting volunteers only, until clinical approval, patient consent, a data-processing agreement and validation on real phones are in place.

The camera can

- Estimate 2D body points, joint angles, left/right symmetry and movement quality.
- Do this only when the needed joints are visible and the set-up checks pass.
- Count repetitions and compare with the target range.

The camera cannot

- Identify the tissue or condition causing pain.
- See bones or vertebrae. That needs X-ray, CT, MRI or ultrasound, read by qualified clinicians.
- Measure real distances (centimetres) from one camera.
- Give disease probabilities, severity grades or diagnoses.

Rules the app always follows

- 1 No number without a good capture.** If capture fails, the value is withheld, the reason is shown and a recapture is requested.
- 2 Sources are never mixed.** Patient reports, camera estimates, algorithmic observations and clinician findings are always labelled and kept apart.
- 3 Every observation shows its working:** measurement, view, confidence, calculation, version and limitations, in the "Why?" view.
- 4 The physiotherapist decides.** Drafts and observations are proposals until a clinician reviews them. Rule sets marked DRAFT still need clinical-lead approval.
- 5 Demo is always labelled.** Simulated data shows a purple "Demonstration data" banner and SIMULATED tags.

For patients

- Exercise in a clear space, with a chair or wall nearby for balance.
- Stop if pain is sharp or getting worse, or if you feel dizzy or unwell.
- The check-in is not monitored in real time. For urgent problems, call your clinic or emergency services.

Troubleshooting

It keeps saying "move closer" or "step back"

Read the Distance line in the checklist: it shows the percentage of the picture your body fills and the range to aim for. Stand so the whole body, head to feet, is in the picture, at the distance the test states (usually 2–3 m). Keep the phone upright. If it still happens, take a screenshot of the checklist and send it to the clinic.

The camera does not start

Allow camera access when the browser asks. If you refused before, allow it in the browser's site settings. Check you gave camera consent in Profile → Privacy & consent. Close other apps that use the camera.

A value shows "—"

A joint the value needs was hidden, out of the picture or unsteady. See "Why a value shows —" in section 08, fix the set-up and recapture.

"Tests are paused"

Your safety answers need your physiotherapist's review first. They will contact you.

The first scan is slow

The body-tracking model downloads once. The Heavy model (Settings) is ~31 MB; use Lite or Full on phones.

The email link says it is invalid or expired

Links work once and expire. Sign in normally, or use Forgot password to get a new link.

I signed up but I don't see the physiotherapist screens

Your email must be on the clinic's physiotherapist list. Ask the clinic owner to add it, then sign out and in again.

"Changes are not yet saved to the clinic server"

You are offline or the server was unreachable. Your changes are kept on the device and sent automatically. A physiotherapist can press Settings → Clinic server → Sync now.

A patient's data is missing on the physio laptop

Check both devices are signed in to real accounts, not demo. Demo data never leaves the device. Wait about 20 seconds or press Sync now.

The report went back to "preliminary"

A capture or review changed after approval. Check the changes and approve again.

Text is too small

Profile → Large text, or use the browser's zoom.

Quick reference cards

Patient: 10 steps

1. Open the site → **Get started** → confirm your email.

2. Setup: About you → consent → your concern → choose an area.

3. **Symptom map**: tap areas, choose symptoms, draw spreading pain.

4. **History**: pain now and worst in 24 h, feel, duration, triggers.

5. **Safety check**: answer honestly.

6. Phone upright on a stand where the test says (usually 2–3 m), whole body in view, light in front.

7. **Movement tests**: wait for all green ticks; recapture any X.

8. **Submit** to your physiotherapist.

9. When your program appears: Home → **Start session**; check in daily.

10. Stop if pain is sharp or worsening. Emergency → call for help.

Physiotherapist: 12 steps

1. Get your email added to the clinic list → sign up or sign in.

2. **Overview**: clear critical alerts and the exception queue first.

3. **Assessments**: open the oldest submitted assessment.

4. History & safety → check flags and the original answers.

5. Captures & replay → Accept / Reject / Repeat each estimate.

6. Left / right or baseline → current.

7. AI draft → Accept / Edit / Reject / Defer each statement.

8. Evidence & reasoning → examine, use "Why?", sign the impression.

9. Report → Approve as clinician-reviewed → download PDF.

10. **Programs** → build, set pain-pause and reassessment → publish.

11. **Schedule** → book visits, mark attended; **Billing** → record payments.

12. End of care → patient → **Discharge** → write, sign, download.

Glossary

TERM	MEANING
Body points (landmarks)	33 points the camera model finds, such as shoulders, hips, knees and ankles.
ROM	Range of motion: how far a joint moves, in degrees.
Camera estimate	A value calculated from body points. It is not a clinical measurement until reviewed.
Algorithmic observation	A camera estimate that crossed a clinic-set threshold. A prompt for review, not a finding.
Valid capture	A recording that passed the quality checks (joints visible, steady, enough frames).
Plumb line	The vertical reference line from the ankles on the posture scan.
FPPA	Frontal plane projection angle: knee alignment seen from the front during a squat.
Clinical impression	The physiotherapist's signed conclusion. Only this appears as a conclusion in the report.
Treatment course	The number of sessions and the fee agreed with a patient; attended sessions are counted from the schedule.
PT code	The patient's short ID in the directory, for example PT-3F9A21.

Dheepika Lab · physiovision-ai-eta.vercel.app · Guide version 6 October 2026. Screens show simulated demonstration data. This guide describes a pilot; it is not clinical advice.